

Anti-pressure Ulcer

TS-10C Alternating Air Mattress & Pump System



(Non-include hospital bed)

Long-term mobility impaired, long-term bedridden, spinal cord injured, disabled after stroke, paralyzed, unconscious in coma, incontinent, diabetic, patients or any patient who has delicate skin, weak physical strength or malnutrition can use the alternating air mattress & pump system to prevent and manage decubitus.



Digital Control Power Unit



Dual External Fuse




Safety Power Management




Wide Elastic Band for Mounting



CPR Quick Deflation Valve

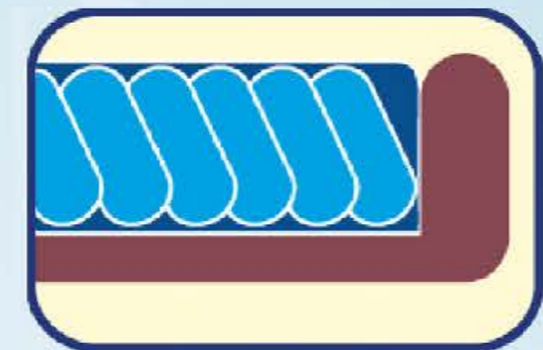

Air Cells Overlap Closely

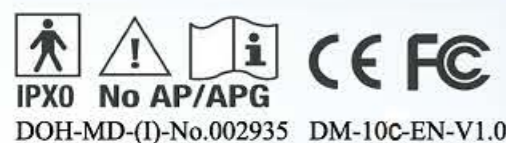

Features of Control Unit: Safety Practical Aesthetic

Internal Pressure Sensor and Feedback System to Monitor the Pressure Variation 24-hours in Real Time



Comfort Setting /Knob Control:

Automatically provides a suitable alternating pressure according to the selected comfort setting. The air mattress can be easily adjusted to a proper hardness/softness ranging from 1 to 10.



Auto Firm Mode (Nursing) Button:

During nursing, it can be difficult to turn a patient and nurse if the mattress is too soft. This can harden the mattress in the quickest way for easy nursing. However, after 30 minutes, the unit will forcefully resume the setting before the button is pressed in order to protect the patient from the risk of developing pressure ulcers because the mattress becomes too hard. Press the button when necessary.



Audible Mute Button:

The alarm becomes audible when Low Pressure Alarm or Power Failure Alarm alerts problem occurrence. The alarm can be muted with this button when it is necessary to turn off the function of alarm sound. To restore the audible alarm function, you have to reset the power switch.



Pressure Monitoring LED Indicator:

When the LED is flashing, this indicates that the pump is in the process of giving sufficient air to the mattress. When the LED is lit on, this indicates that the pump is inflating the mattress according to the selected comfort setting. If the LED is not lit on, this indicates that the mattress has reached to an air pressure according to the selected comfort setting.

- Uniquely and exclusively designed high quality Knob are in good taste and aesthetic. The texture and protruded positioning points at the sides allow the air mattress to be easily adjusted to a proper softness/hardness from 1 to 10 number.
- 3-PIN separable power plug design can avoid the electrical hazard of the power cord while it is being pulled out during emergency.
- Safe external dual-fuse design enhances the protection ratings and allows the fuses to be replaced directly without the inconvenience of disassembly of the machine.
- Patented internal mechanism is featured with low vibration, low noise, low power consumption, and modular design for convenient production and maintenance.
- Entire unit is designed with streamline outline, lightweight body, perfect geometric ratio, and visual comfort appearance which provide dignified care for the user.
- The compact size feature: The system is of a compact size, making it easy for placement and storage and taking into account the streamlined system appearance and practicality.
- Dual-hook design reinforces the hanging strength and can be adjusted at different hanging angles. Its rigid steel components can greatly endure the external pulling forces.
- Unit can be conveniently operated by being placed vertically, horizontally or hung on a hook.
- Filter sponge can be replaced externally without the trouble of disassembly for maintenance.

ABC Alternating



Features of Air Mattress: Comfort Ease Durability

- 4" Air cells – 20pcs totally :Entire mattress is constructed with a large number of air cells and a denser spacing between air cells that the body weight can be supported more evenly and more comfortable.
- Head End 3 Hoses No Alternating Design.
- External connecting tube are protected with dust-proof sleeve so as to avoid the accumulation of dusts and stains on the surface.
- Patent Design of CPR Quick Valve are equipped for fast deflation of the air mattress so as to facilitate the cardiopulmonary resuscitation (CPR) during the emergency.
- Rear side of the air mattress is equipped with wide elastic bands on the four corners for simple and easy mounting on the original foam mattress or at the frame of the medical bed.

MODEL: TS-10C

Control Unit Warranty:1 year	Input Power:	AC100-120V/60Hz or AC220-240V/50Hz
	Max. Power Consumption:	10W max
	External Fuse:	T1A1/250V
	Dimensions:	L26*W8*H13 cm
	Weight:	Approx. 1.5 kg
	Air Flow Rate:	Approx. 8LPM
	Alternating Mode:	ABC Alternating
Air Mattress Warranty:1 year (air cells only)	Cycle time:	10Minutes
	Appliance Category:	Class I, Type BF
	Certification:	CE , EN60601,EN60601-1(MDD), EN 60601-1-2 , FCC Part 18-Class B
	Dimensions:	L190*W85*H10 cm
	Weight:	Approx. 3.8 kg
	Load:	30-120 kg
	Air cells Materials:	Nylon + PVC
Air Mattress	Base Materials:	Nylon + PVC
	Mattress Cover :	Nylon/PU Coating
	Made in Taiwan	*The product specifications may be subject to change according to the market without further notice.



Manufacturer: TRUE SOURCE Technology Co., Ltd.
4F, No.3, Ln. 130, Minquan Rd., Xindian Dist., New Taipei City 231, Taiwan (R.O.C.)
TEL: +886-2-86671780 FAX: +886-2-86671782
<http://www.medical-truesource.com>

Distributor

Dignified Healthcare  Thoughtful Concern